


RX Cytology Brush

Wireguided Cytology Brush

Boston
Scientific

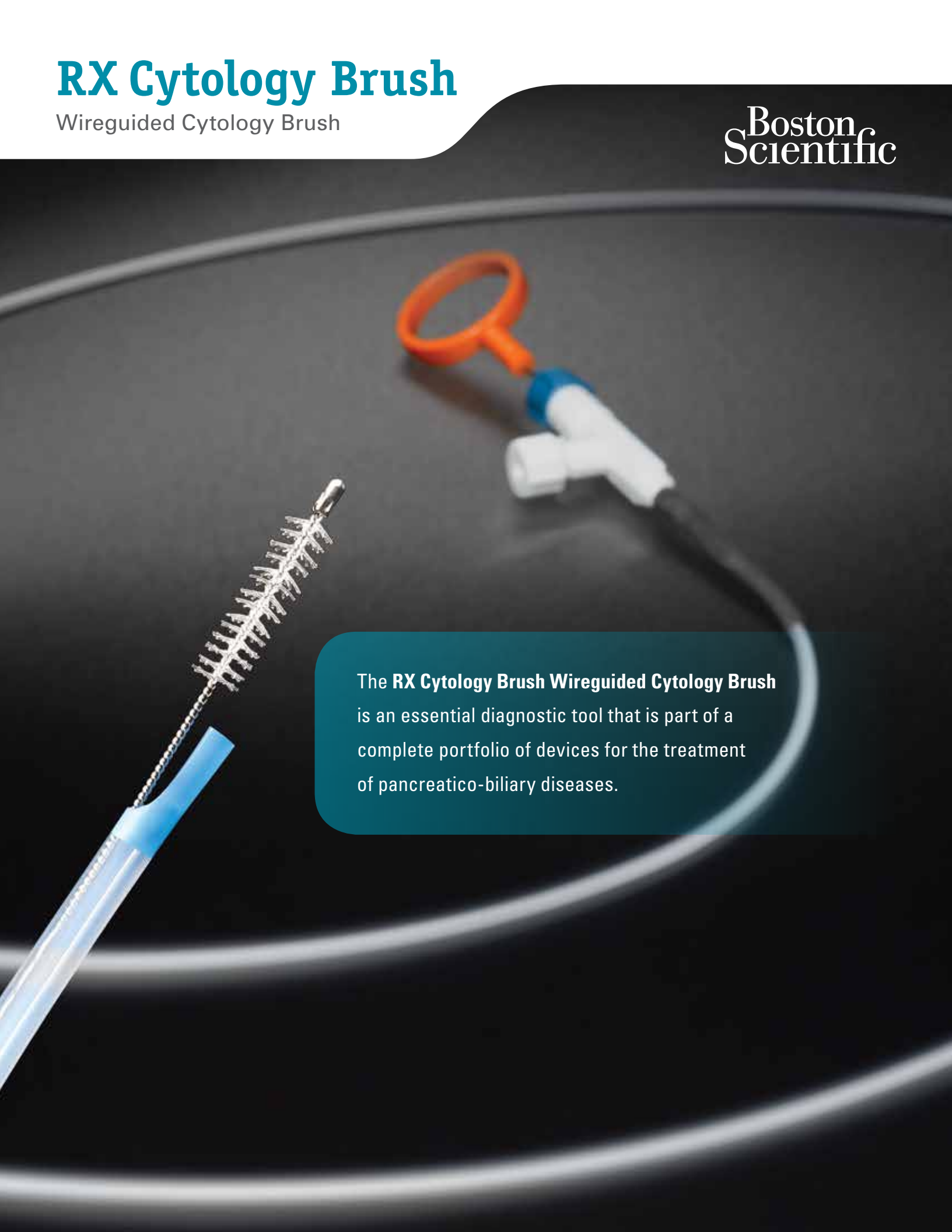

The **RX Cytology Brush Wireguided Cytology Brush** is an essential diagnostic tool that is part of a complete portfolio of devices for the treatment of pancreatobiliary diseases.

RX Cytology Brush


Wireguided Cytology Brush



The RX Cytology Brush Wireguided Cytology Brush is a double-lumen wireguided catheter for the collection of cells and the detection of malignant neoplasms in the biliary duct.

DEVICE FEATURES

- Handle has reinforced metal tube for reliable actuation.
- Distal radiopaque marker aids in catheter visualization under fluoroscopy for proper positioning.
- Open channel allows for quick exchange of the guidewire.

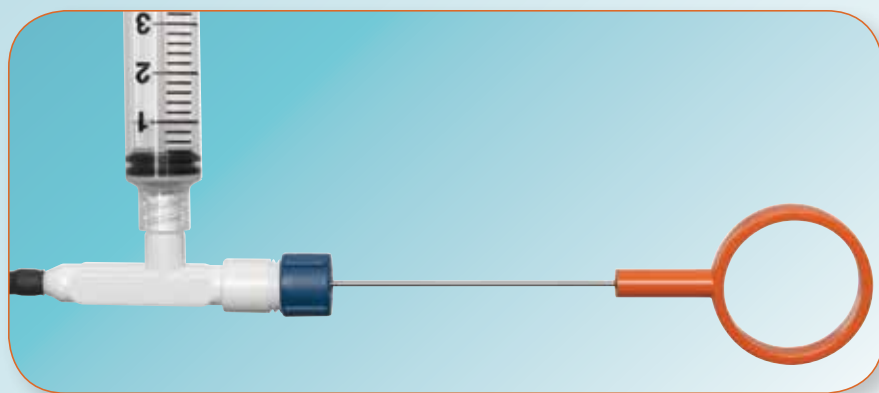

Brush lumen may be used to inject contrast dye through the proximal luer.

ORDERING INFORMATION

Order Number	Product Description	Brush Width (mm)	Catheter O.D. (Fr)	Working Length (cm)	Recommended Guidewire (in)	Units
M00545000	RX Cytology Brush	2.1	8	200	.035	Each

Indications, Contraindications, Warnings and Instructions for Use can be found in the product labeling supplied with each device.

Caution: Federal (U.S.) law restricts this device to sale by or on the order of a physician.

CAUTION: The law restricts these devices to sale by or on the order of a physician.

**Boston
Scientific**

Advancing science for life™

Boston Scientific Corporation
300 Boston Scientific Way
Marlborough, MA 01752-1234
www.bostonscientific.com/endo-resources

Ordering Information
1.888.272.1001

©2014 Boston Scientific Corporation
or its affiliates. All rights reserved.

ENDO-253201-AA July 2014